
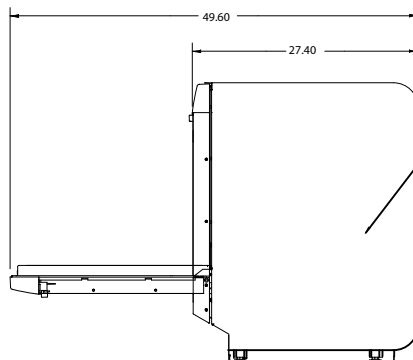


Electrical Requirements:  
115 volts, 60 hz, 16  
amps

Interior

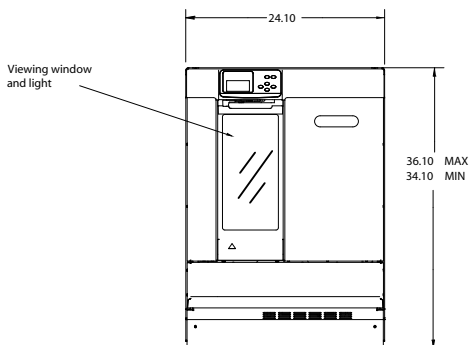

Exposed insulation  
backing  
Typical on all  
undercounter models

Exterior

## GENERAL DESCRIPTION

### 4420330 UNDERCOUNTER FLASKSCRUBBER WASHER

#### GENERAL FEATURES

- Lower spindle rack, of type 304 stainless steel, with 36 detachable large spindles for washing narrow glassware. Eighteen glassware holders, 36 clips and 10 additional interchangeable small spindles are included. Rack accommodates a variety of inserts for holding breakers, culture tubes and other specialty glassware (inserts sold separately).
- Water temperatures reach 82° C (180° F) on 115v models to sanitize glassware and for enhanced washing and faster drying.
- LCD information center for customizing cycles and monitoring operation, including internal water temperature (°C or °F).
- Water recalculation rates of 96 gal/min (363 liters) when operated on 115 volt, 60hz for thorough cleaning.
- Built in purified water pump to bring non-pressurized or pressurized purified water into the tank for up to six pure water rinses.
- Built in steam generator that produces hot vapor to penetrate and remove dried residue for cleaner glassware.
- Manual detergent dispenser for powder or liquid detergent.
- Manual rinse aid solution dispenser.
- Built-in forced air drying up to 99 minutes and programmable from 38-70°C (100-158°F).
- Factory cycle programs. All 115 volt models have seven programs: Rinse only, Plastic, Glass, Glass Plus, Scientific, Scientific Plus, and Dry only. All 230 v models have three additional programs for a total of ten: Intense, Intense Plus, and Extreme.
- Alarms with display for HOT GLASS, PLEASE WAIT, WATER HIGH, and OVERFLOW.
- Two user cycle programs.
- User delayed start (up to 8 hours).
- Type 304, brushed stainless steel door, tank upper and lower wash arms. No exterior sides or top.
- One piece, stainless steel, fine mesh particle drain filter to trap debris in tank to protect the pump.
- No center water distribution tower allowing for greater glassware capacity.
- Four leveling feet.
- 3/8" female NPT inlet fitting for hot tap water connection.
- Factory installed drain hose and band clamp.
- Two year warranty on parts and labor.

#### requirements

- Minimum inlet water temperature 49°C (120°F)
- Recommended inlet water temperature 60° (140°F)
- Minimum purified water pressure 0 psi
- Water consumption per fill 3.4 gallons (13 liters)
- Hot tap water pressure at the washer 20-120 psi
- Hot tap water flow rate 1.25 gallons (4.7 liters) per minute.
- 3/8" OD copper tubing for connection to the water inlet valve.
- Purified water with supply piping minimum 3/8" ID to permit at least 0.9 gallon (3.4 liters) per minute flow.

Glassware Washer  
for:



Labconco Corporation  
10000 E. 15th Avenue  
Kansas City, MO 64132-2496  
Phone: 816-254-2000  
www.labconco.com

No.	Date

Date:  
09/26/07

Drawing #  
4420330

Drawn By:

Checked By:

Approved By:

Sheet No.

Project No.